



National Institute of Standards & Technology

Certificate

Standard Reference Material[®] 2241

Relative Intensity Correction Standard for Raman Spectroscopy: 785 nm Excitation

This Standard Reference Material (SRM) is a certified spectroscopic standard for the correction of the relative intensity of Raman spectra obtained with instruments employing 785 nm laser excitation. A unit of SRM 2241 consists of an optical glass that emits a broadband luminescence spectrum when excited at this laser wavelength. The measurements of the relative spectral intensity are traceable to a lamp-illuminated integrating sphere source which has been calibrated for its spectral radiance at NIST. The shape of the mean luminescence spectrum of this glass is described by a mathematical expression that relates the relative spectral intensity to the wavenumber (cm^{-1}) expressed as the Raman shift from the excitation laser wavelength. This model, together with a measurement of the luminescence spectrum of the standard, can be used to determine the spectral intensity-response correction that is unique to each Raman system. The resulting instrument-intensity-response correction may then be used to obtain Raman spectra that are largely free from instrument-induced spectral artifacts.

Certified Values: The polynomials describing the relative luminescence spectrum of SRM 2241 and associated confidence curves are given in Table 1. The measurand is relative luminescence measured as a function of Raman shift (cm^{-1}) from the laser excitation wavelength of 785 nm. Metrological traceability is to the NIST spectral radiance scale. A NIST certified value is a value for which NIST has the highest confidence in its accuracy in that all known or suspected sources of bias have been investigated or taken into account [1].

Expiration of Certification: The certification of **SRM 2241** is valid, within the measurement uncertainty specified, until **30 November 2020**, provided the SRM is handled and stored in accordance with the instructions given in this certificate (see “Instructions for Handling, Storage, and Use).” The certification is nullified if the SRM is damaged, contaminated, or otherwise modified.

Maintenance of SRM Certification: NIST will monitor this SRM over the period of its certification. If substantive technical changes occur that affect the certification before the expiration of this certificate, NIST will notify the purchaser. Registration (see attached sheet or register online) will facilitate notification.

Production and certification of this SRM were performed by S.J. Choquette and A.A. Urbas of the NIST Biosystems and Biomaterials Division. The SRM units were cut and polished by J. Fuller, formerly of NIST.

Statistical consultation was provided by S.D. Leigh formerly of the NIST Statistical Engineering Division.

Support aspects involved in the issuance of this SRM were coordinated through the NIST Office of Reference Materials.

Preparation and certification of this SRM were supported in part by the Test & Evaluation and Standards Division, Science and Technology Directorate, of the Department of Homeland Security.

Steven J. Choquette, Acting Chief
Biosystems and Biomaterials Division

Gaithersburg, MD 20899
Certificate Issue Date: 24 August 2015
Certificate Revision History on Last Page

Robert L. Watters, Jr., Director
Office of Reference Materials

Certified Values: The certified values of the coefficients of the model describing the mean shape of the luminescence spectrum of SRM 2241, excited at 785 nm, are listed in Table 1. The spectrum and its associated expanded uncertainty [$\pm 95\%$ confidence curves (CC)] are shown in Figure 1 [2,3]. This model is certified to describe the luminescent response of the SRM when it is measured in the temperature range of 20 °C to 25 °C. This model certifies the shape of the luminescence spectrum between 200 cm^{-1} and 3500 cm^{-1} Raman shift for excitation with a 785 nm laser.

Certification Uncertainty: The coefficients of the models that express the expanded uncertainty ($\pm 95\%$ confidence curves) of the certified model of $I_{\text{SRM}}(\Delta\nu)$ are listed in Table 1. These models are used to calculate the upper and lower 95 % confidence curves for the mean luminescence spectrum of SRM 2241. These upper and lower bounds are shown in Figure 1.

The confidence curves were calculated by applying the point-by-point uncertainty analysis of [2,3] to the data from which the certified model was derived. The curves represent the uncertainties on the mean luminescence spectrum of SRM 2241 based on measurements using three different spectrometers. Components of the uncertainty include the uncertainty in the calibration source and various estimable systematic errors from the operation of the Raman instruments used in the measurements. Careful measurements of the glass have shown it to be spatially homogeneous in spectral luminescence. No significant changes in the shape of the luminescence spectrum occur over the range of laser power densities commonly used in Raman instruments.

Table 1. Coefficients of the Certified Polynomial ^(a) and of the 95 % Confidence Curves^(b)
Describing the Mean Luminescence Spectrum of SRM 2241 for 785 nm excitation.
(Valid for Temperatures of 20 °C to 25 °C)

Polynomial Coefficient	Certified Value Polynomial Coefficient ^(c) 20 °C to 25 °C	Polynomial Coefficient ^(c) of the $\pm 95\%$ Confidence Curves ^(b)	
		+95 % CC	−95 % CC
A ₀	9.71 937 E−02	1.04 276 E−01	9.01 111 E−02
A ₁	2.28 325 E−04	2.39 131 E−04	2.17 519 E−04
A ₂	−5.86 762 E−08	−7.81 489 E−08	−3.92 035 E−08
A ₃	2.16 023 E−10	2.32 243 E−10	1.99 803 E−10
A ₄	−9.77 171 E−14	−1.03 769 E−13	−9.16 653 E−14
A ₅	1.15 596 E−17	1.23 774 E−17	1.07 417 E−17

^(a) A NIST certified value represents a value for which NIST has the highest confidence in its accuracy in that all known or suspected sources of bias have been fully investigated or taken into account [1].

^(b) The consensus curve is a point-by-point weighted mean of the average responses of three instruments [4,5], fitted by the polynomial model. The uncertainty curves are polynomial models of point-by-point expanded uncertainties, with coverage factor $k = t_{0.975,2} = 4.303$ (95 % confidence), calculated by combining a between-instrument variance incorporating interinstrumental Type B uncertainties with a pooled within-method variance [3], following the ISO/JCGM Guide [2].

^(c) Where $I_{\text{SRM}}(\Delta\nu) = A_0 + A_1 \times (\Delta\nu)^1 + A_2 \times (\Delta\nu)^2 + A_3 \times (\Delta\nu)^3 + A_4 \times (\Delta\nu)^4 + A_5 \times (\Delta\nu)^5$,
for $\Delta\nu = 200 \text{ cm}^{-1}$ to 3500 cm^{-1} Raman Shift relative to 785 nm excitation.
 $I_{\text{SRM}}(\Delta\nu)$ has units of photons $\text{s}^{-1} (\text{cm}^{-2}) (\text{cm}^{-1})^{-1}$

Physical Description: SRM 2241 is a chromium-doped (mole fraction of 0.02 % Cr_2O_3) sodium borosilicate matrix glass. Each unit of this SRM consists of a glass slide that is approximately 10 mm in width \times 10 mm in length \times 1.65 mm in thickness, with one surface optically polished and the opposite surface ground to a frosted finish using a 400 grit polish. Two mounts are furnished with the slide. One is a 12.5 mm square cuvette-style optical mount. This mount is designed for the typical 12.5 mm sampling accessories widely used in chemical spectroscopy, i.e., absorbance, fluorescence, etc. The glass slide is retained, frosted side out, in a slot on the front face of the holder. Two plastic springs, also retained in the slot, hold the glass slide in place while allowing for positioning within the slot to accommodate different beam heights. Grooves on the sides of the slot require that the glass slide and springs be loaded into the slot from the bottom of the holder. The other mount is a 2.5 cm \times 7.6 cm \times 0.3 cm microscope slide-style holder. On the top face is a rectangular slot to retain the glass slide over a circular aperture in the center of the holder. Two plastic springs serve to hold the glass slide in the slot centered over the circular aperture. Removal of the SRM glass for measurements that are physically hindered by the holders do not alter the certified properties of this SRM.

Measurement Conditions: The certification measurements of the luminescence spectrum of SRM 2241 were made using three spectrometer systems: one commercial Raman microscopy system operated in a 180° backscatter geometry, one Raman system based on 460 mm focal length spectrograph designed for array detectors operated in a 180° backscatter geometry, and one commercial near-IR spectrofluorometer operated in a 90° geometry [6]. Excitation laser wavelength was measured daily using a wavemeter. Excitation laser sources included an argon ion laser pumped titanium:sapphire solid-state tunable laser (785.00 nm \pm 0.01 nm) and a frequency-stabilized 785 nm diode laser source (784.98 nm \pm 0.02 nm). Either edge or notch Rayleigh rejection filters were used on all systems. The absolute wavenumber or wavelength axis of each spectrometer was calibrated using emission lines from low-pressure pen lamps operated with a DC power supply. The y-axis (relative spectral intensity) of each system was calibrated with a white-light, uniform-source, integrating sphere that had been calibrated for irradiance at NIST or by a high-emissivity blackbody radiator with furnace temperature calibrated by the manufacturer. All certification data were acquired at nominal room temperature (22 °C).

INSTRUCTIONS FOR HANDLING, STORAGE, AND USE

Handling and Storage: The glass slide should be handled only in its mount. Although not recommended, the glass standard may be removed from its mount without altering the certified properties of the glass. When not in use, the SRM should be stored in the container provided or in one providing comparable mechanical protection.

Use: SRM 2241 is used to provide Raman spectra corrected for relative intensity and is not intended for use as a standard for the determination of absolute spectral irradiance or radiance. This requires a measurement of its luminescence spectrum on the Raman instrument and then a mathematical treatment of both this observed luminescence spectrum, and the observed Raman spectrum of the sample. This SRM is intended for use in measurements over the temperature range of 20 °C to 25 °C. It may also be used for Raman excitation with lasers that range from 784 nm to 786 nm in excitation wavelength. Use of this SRM at temperatures outside of the certification temperature range is not currently supported.

For proper use of this procedure, attention must be paid to the following experimental conditions. The spectrometer laser and x-axis should be calibrated using the manufacturer's recommended methods. Validation of the Raman shift axis may be accomplished by referring to ASTM E1840-96 [7]. It is important that the laser excitation be incident only on the frosted surface of the glass. The shape of the spectral luminescence will have some sensitivity to the placement of the glass surface relative to the collection optics of the spectrometer, which is minimized by scattering from the frosted surface. Measurement conditions should be arranged to furnish a spectrum of optimum signal-to-noise ratio of the SRM. The luminescence spectrum must be acquired over the same Raman range as that of the sample.

The relative intensity of the measured Raman spectrum of the sample can be corrected for the instrument-specific response by a computational procedure that uses a correction curve. This curve is generated using the certified model and the measured luminescence spectrum of the SRM glass. For the spectral range of certification, $\Delta\nu = 200 \text{ cm}^{-1}$ to 3500 cm^{-1} , compute the elements of the certified relative mean spectral intensity of SRM 2241, $I_{\text{SRM}}(\Delta\nu)$, according to equation (1)

$$I_{\text{SRM}}(\Delta\nu) = A_0 + A_1 \times (\Delta\nu)^1 + A_2 \times (\Delta\nu)^2 + A_3 \times (\Delta\nu)^3 + A_4 \times (\Delta\nu)^4 + A_5 \times (\Delta\nu)^5 \quad (1)$$

where $(\Delta\nu)$ is the wavenumber in units of Raman shift (cm^{-1}) and the A_n 's are the coefficients listed in Table 1. The elements of $I_{\text{SRM}}(\Delta\nu)$ are obtained by evaluating equation (1) at the same data point spacing used for the acquisition of the luminescence spectrum of the SRM and of the Raman spectrum of the sample. $I_{\text{SRM}}(\Delta\nu)$ has been normalized to unity and is a relative unit expressed in terms of photons $\text{s}^{-1} \text{ cm}^{-2} (\text{cm}^{-1})^{-1}$. The data sets that are the measured glass luminescence spectrum, S_{SRM} , and the measured Raman spectrum of the sample, S_{MEAS} , must have the units of Raman shift (cm^{-1}). The elements of the correction curve $I_{\text{CORR}}(\Delta\nu)$, defined by equation (2), are obtained from $I_{\text{SRM}}(\Delta\nu)$ and the elements of the glass luminescence spectrum, $S_{\text{SRM}}(\Delta\nu)$, by

$$I_{\text{CORR}}(\Delta\nu) = I_{\text{SRM}}(\Delta\nu) / S_{\text{SRM}}(\Delta\nu). \quad (2)$$

The elements of the intensity-corrected Raman spectrum, $S_{\text{CORR}}(\Delta\nu)$, are derived by multiplication of the elements of the measured Raman spectrum of the sample, $S_{\text{MEAS}}(\Delta\nu)$, by the elements of the correction curve [8]

$$S_{\text{CORR}}(\Delta\nu) = S_{\text{MEAS}}(\Delta\nu) \times I_{\text{CORR}}(\Delta\nu). \quad (3)$$

The certified model, equation (1), is certified for use **between 200 cm^{-1} and 3500 cm^{-1}** . The certified model is intended as a simple numerical descriptor of the spectral response observed over the wavenumber range studied. It is not intended as a physically meaningful model. **The model coefficients listed in Table 1 cannot be used to extrapolate the limits of certification without incurring significant error. Extrapolation of the model outside the certification limits of 200 cm^{-1} and 3500 cm^{-1} is not a supported use of this SRM.**

Luminescence Spectrum on the Wavelength Scale: The equation describing the mean luminescence spectrum of the glass SRM is given in equation (1), where $(\Delta\nu)$ is the wavenumber in units of Raman shift (cm^{-1}). For correction of spectra where the x-axis is in wavelength with units of nanometers, the same model coefficients can be used to calculate $I_{\text{SRM}}(\lambda)$ through the following coordinate transformation:

$$I_{\text{SRM}}(\lambda) = [10^7/\lambda^2] \times [A_0 + A_1 \times Z^1 + A_2 \times Z^2 + A_3 \times Z^3 + A_4 \times Z^4 + A_5 \times Z^5] \quad (4)$$

where

$$Z = 10^7 \times [(1.0/\lambda_L) - (1.0/\lambda)] \quad (5)$$

and λ_L is the wavelength of the laser in nanometers and λ is the wavelength in nanometers. The prefactor of 10^7 in the first term of equation (4) is needed only if it is desired to preserve the numerical value of spectral areas computed relative to the two x-axis coordinate systems.

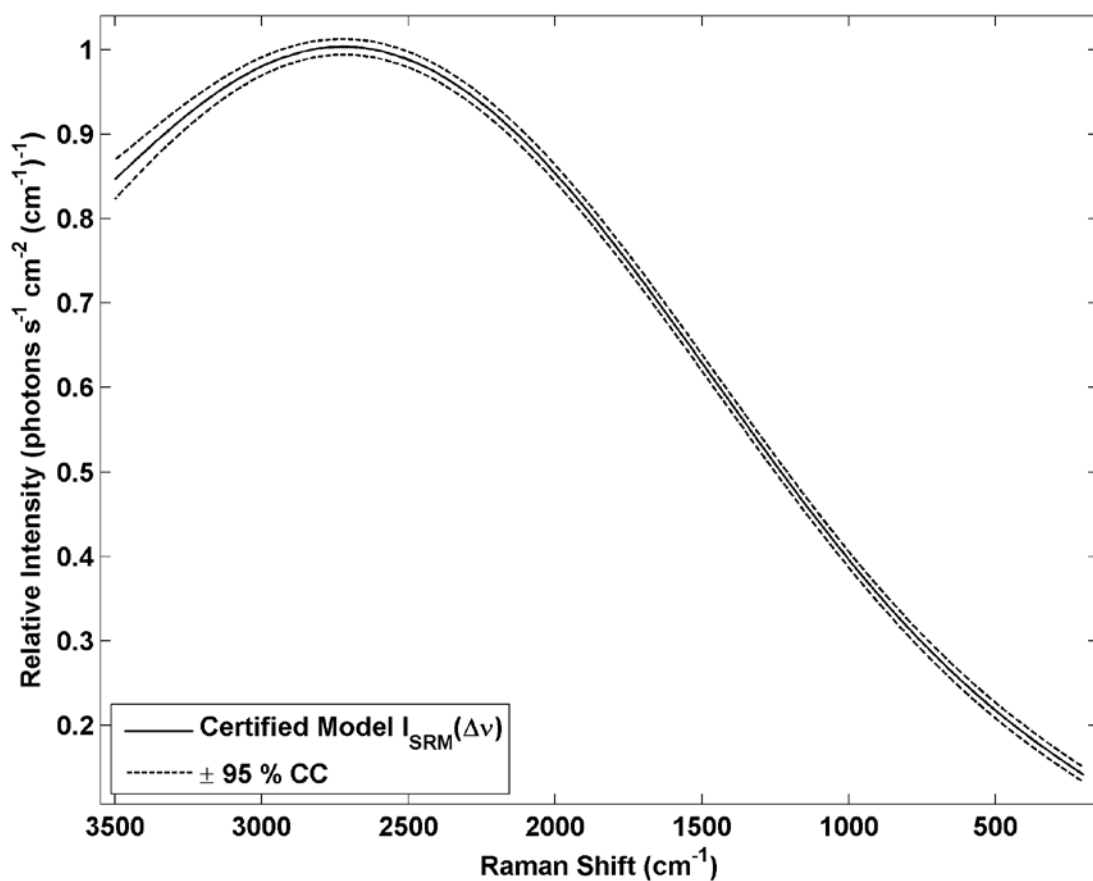


Figure 1. Fifth-order certified polynomial for SRM 2241 when excited at 785 nm. The x-axis is expressed in Raman shift (cm^{-1}). The y-axis is on a relative scale and normalized to unity with the dimensions of photons $s^{-1} cm^{-2} (cm^{-1})^{-1}$. Dashed lines denote 95 % confidence curves.

REFERENCES

- [1] May, W.; Parris, R.; Beck II, C.; Fassett, J.; Greenberg, R.; Guenther, F.; Kramer, G.; Wise, S.; Gills, T.; Colbert, J.; Gettings, R.; MacDonald, B.; *Definition of Terms and Modes Used at NIST for Value-Assignment of Reference Materials for Chemical Measurements*; NIST Special Publication 260-136 (2000); available at <http://www.nist.gov/srm/publications.cfm> (accessed Aug 2015).
- [2] JCGM 100:2008; *Evaluation of Measurement Data — Guide to the Expression of Uncertainty in Measurement* (GUM 1995 with Minor Corrections); Joint Committee for Guides in Metrology (2008); available at http://www.bipm.org/utls/common/documents/jcgm/JCGM_100_2008_E.pdf (accessed Aug 2015); see also Taylor, B.N.; Kuyatt, C.E.; *Guidelines for Evaluating and Expressing the Uncertainty of NIST Measurement Results*; NIST Technical Note 1297; U.S. Government Printing Office: Washington, DC (1994); available at <http://www.nist.gov/pml/pubs/index.cfm> (accessed Aug 2015).
- [3] Horn, S.D.; Horn, R.A.; Duncan, D.B.; *Estimating Heteroscedastic Variances in Linear-Models*; J. Am. Stat. Assoc., Vol. 70(350), pp. 380–385 (1975).
- [4] Rukhin, A.L.; *Weighted Means Statistics in Interlaboratory Studies*; Metrologia, Vol. 46, pp. 323–331 (2009).
- [5] Dersimonian, R.; Laird, N.; *Meta-Analysis in Clinical Trials*; Control. Clin. Trials, Vol. 7(3), pp. 177–188 (1986).
- [6] Choquette, S.J.; Etz, E.S.; Hurst, W.S.; Blackburn, D.H.; Leigh, S.D.; *Relative Intensity Correction of Raman Spectrometers: NIST SRMs 2241 Through 2243 for 785 nm, 532 nm, and 488 nm/514.5 nm Excitation*; Appl. Spectrosc., Vol. 61, pp. 117–129.
- [7] ASTM E1840–96, *Standard Guide for Raman Shift Standards for Spectrometer Calibration*; Annu. Book ASTM Stand., Vol. 03.06 (2007).
- [8] Frost, K.J.; McCreery, R.L.; *Calibration of Raman Spectrometer Response Function with Luminescence Standards: An Update*; Appl. Spectrosc., Vol. 52(12), pp. 1614–1618 (1998).

Certificate Revision History: 24 August 2015 (Change of expiration date; editorial changes); 28 February 2014 (Update of storage information; editorial changes); 02 December 2010 (SRM 2241 Recertification; Extension of certification period; editorial changes.); 30 May 2008 (Extension of certification period); 08 May 2007 (Update of expiration date); 27 June 2002 (Original certificate date).

Users of this SRM should ensure that the Certificate in their possession is current. This can be accomplished by contacting the SRM Program: telephone (301) 975-2200; fax (301) 348-3730; e-mail srminfo@nist.gov; or via the Internet at <http://www.nist.gov/srm>.